



# 2D Heterostructure Preparation




**Jae-hoon Ji**



# Stamp Preparation

## A. Prepare PDMS stamp

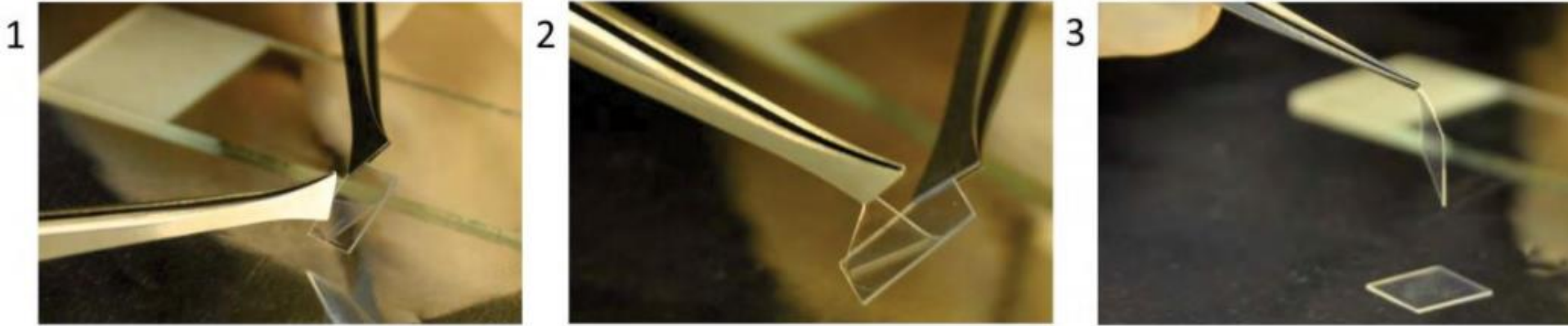


- Material (\$ 120)
  - PDMS material - Sylgard 184 (Dow Corning, USA), \$ 45
  - PPC material - Propylene carbonate (Sigma-Aldrich, CAS 25511-85-7) \$ 75
- Equipment (\$ 306, used for releasing bubbles inside PDMS to reduce surface roughness)
  - Vacuum Desiccator (Southern Labwaer) \$ 56
  - Oil free vacuum pump (Rocker) \$ 250



# Stamp Preparation

## A. Prepare PDMS stamp

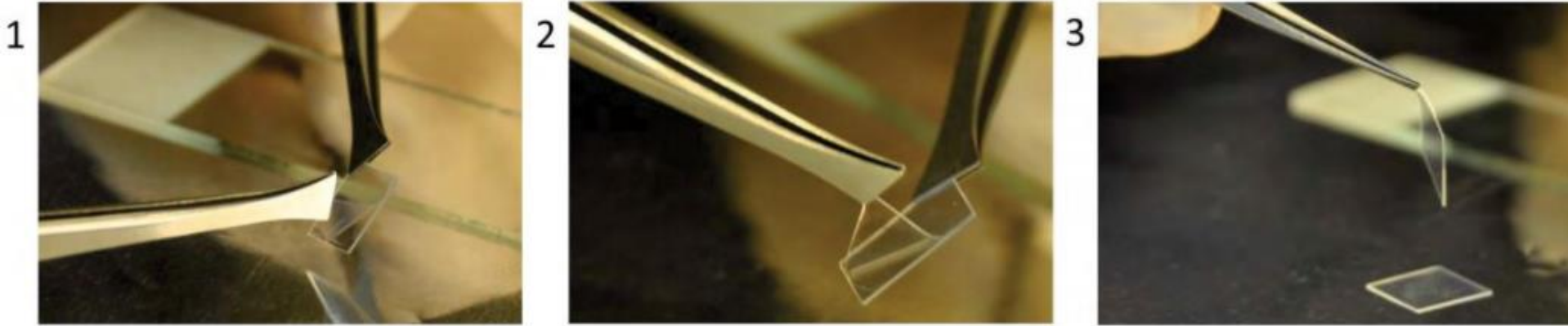

- PDMS preparation


- Prepare the pre-polymer base without bubble (Under vacuum condition for 0.5 hour)
- Mix pre-polymer base and cross linking curing agent (10:1~15:1)
- Pour the mixed solution in the container
- Bubble releasing (Under vacuum condition for 1 hour)
- Cure the solution for 4 hour, 4 C
- If the surface of the stamp is not clean enough => Plasma treating

- PPC preparation


- Coat the surface of the stamp

# Stamp Preparation

## B. Attach the stamp on the glass

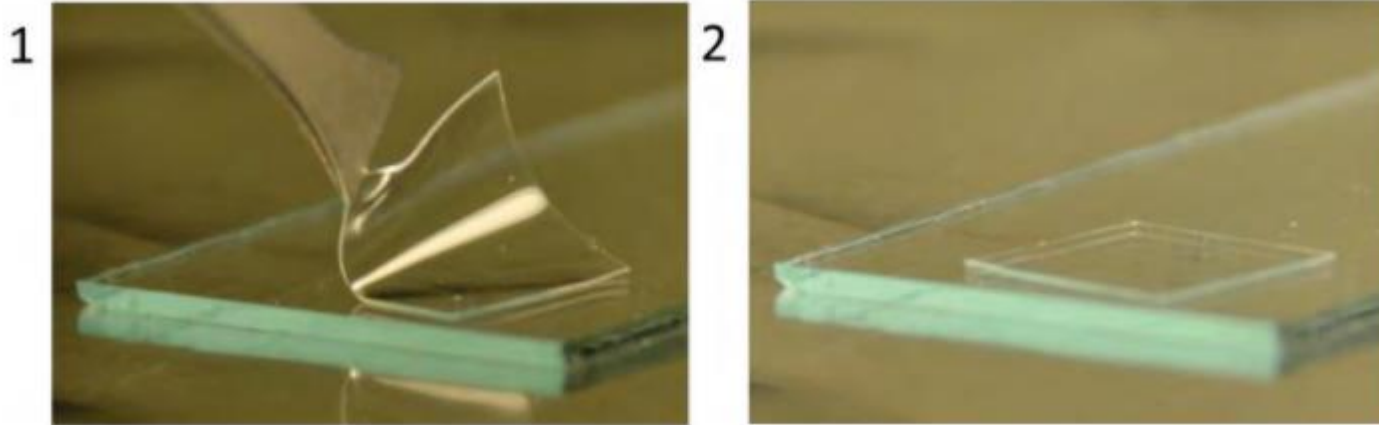



## C. Transfer TMDC flakes on the stamp

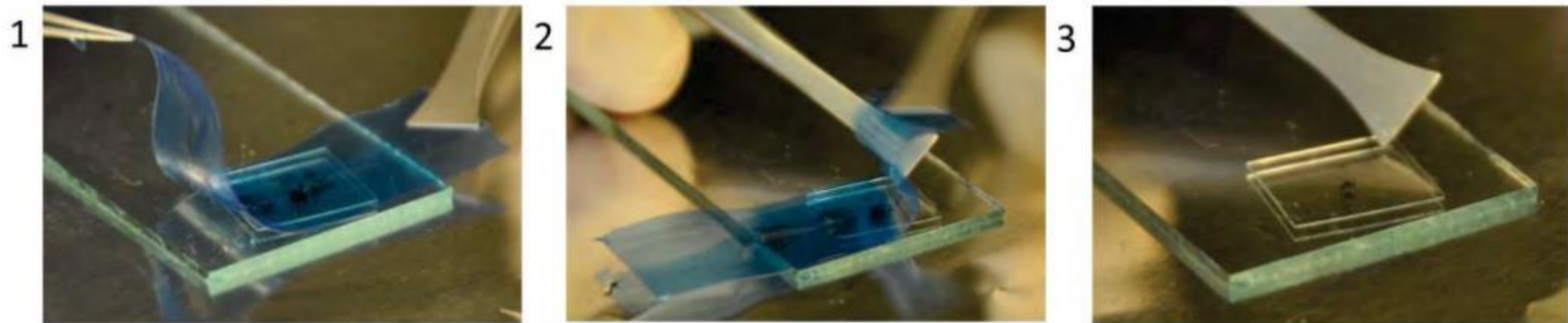



# Transfer method using micromanipulator

## D. Transfer TMDC flakes on the substrate using micromanipulator and optical camera

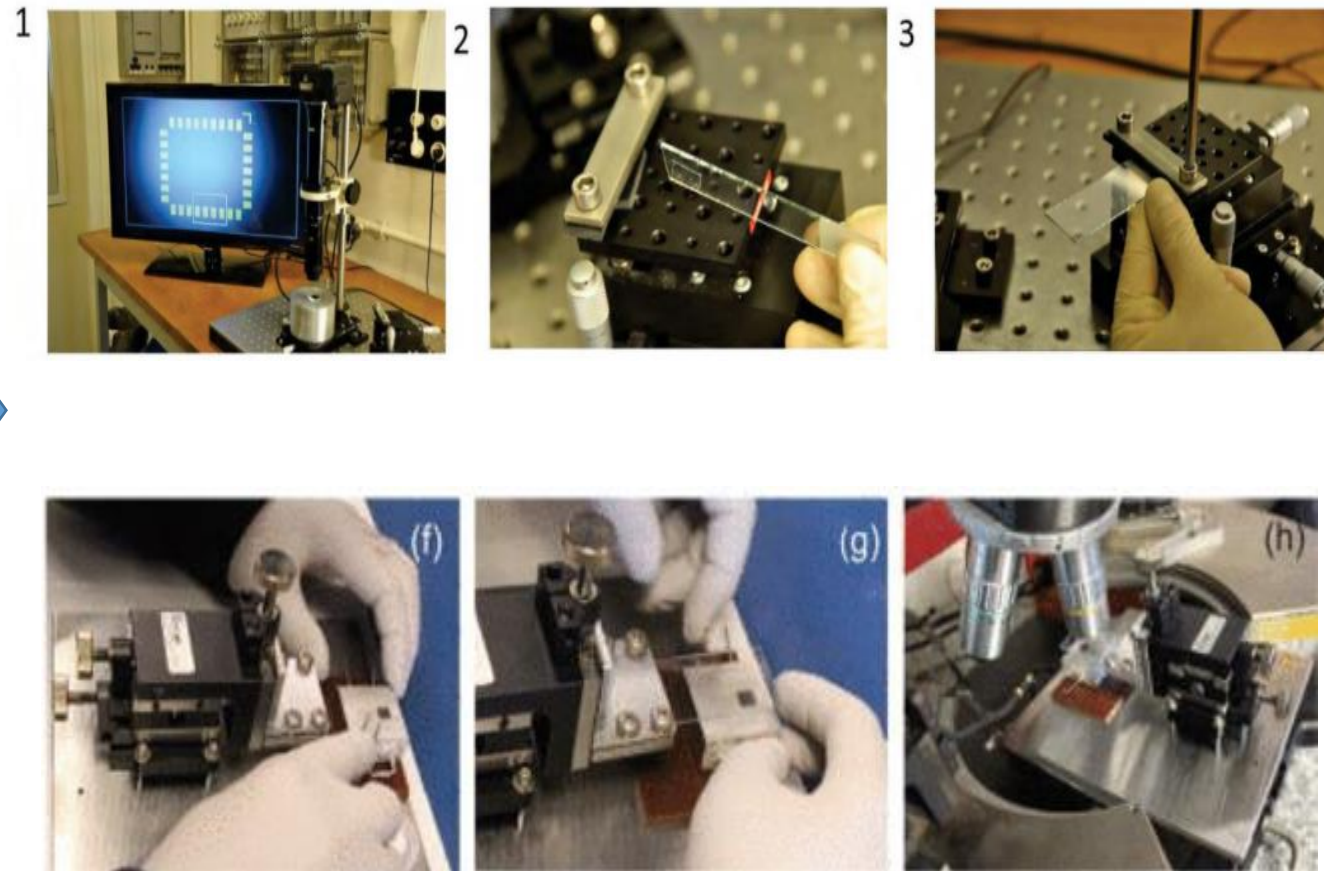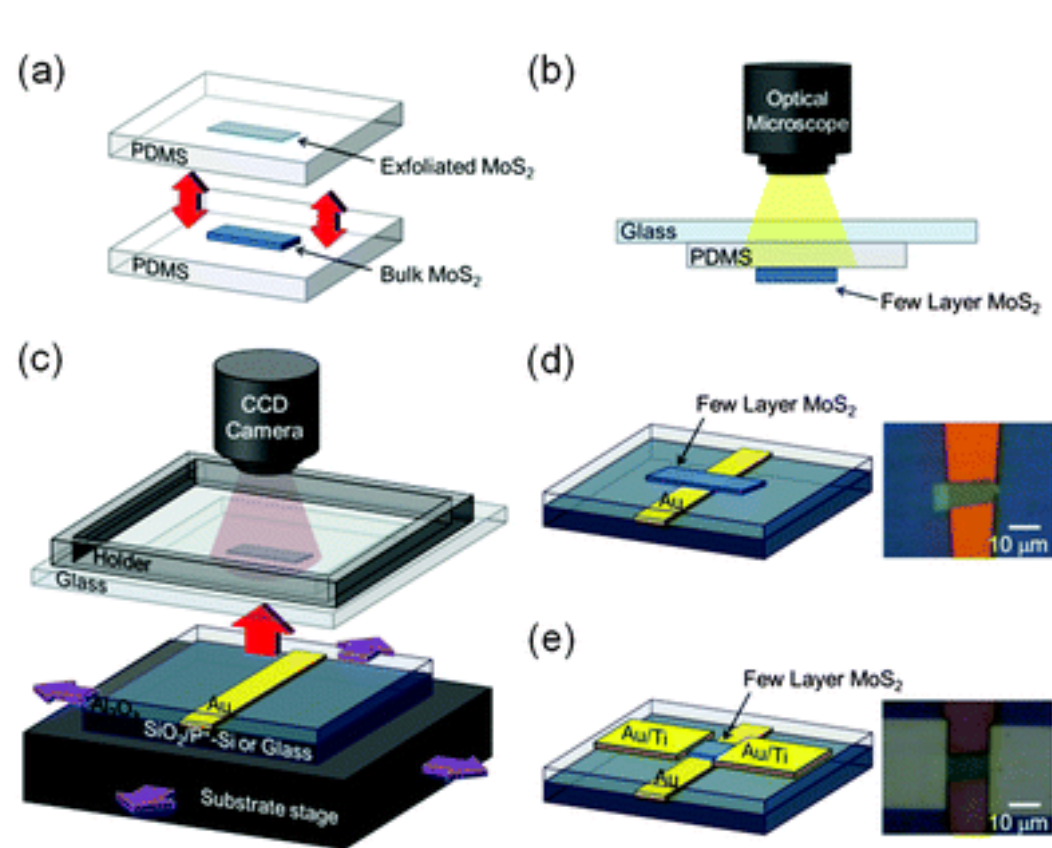

# Transfer Equipment


- Optical scope (AMScope, \$850)



## Micropositioner



MM1-3

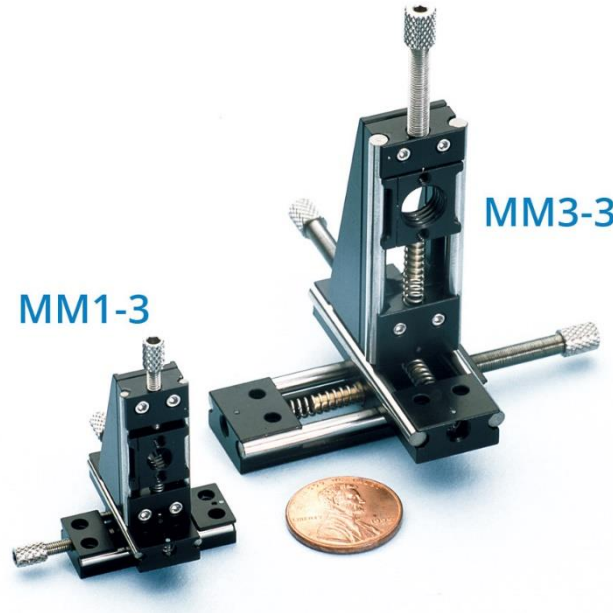


- up to 225x magnification
- Have to find fairly large flakes
- Ask the provider whether we can get more high magnification range.

# Transfer Equipment

## (a) Using Long Working Distance Objective

(a)

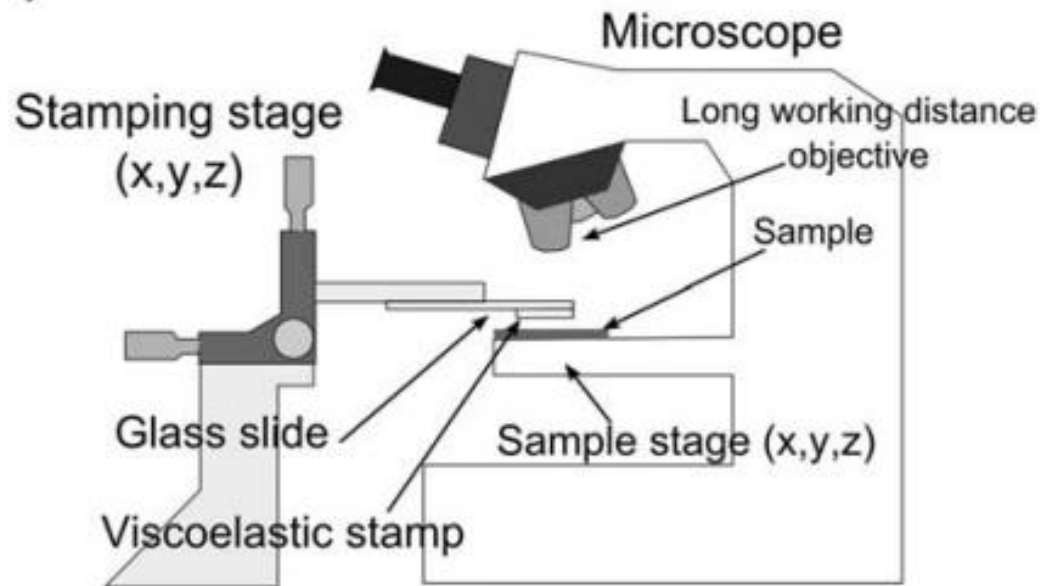

Objective – too expensive

+1	Qty	Docs	Part Number - Price	Available / Ships
	<input type="text"/>		<a href="#">MY5X-802</a> \$695.00	✓ 3-5 Days
	<input type="text"/>		<a href="#">MY10X-803</a> \$873.00	✓ Today
	<input type="text"/>		<a href="#">MY20X-804</a> \$2,056.00	Lead Time
	<input type="text"/>		<a href="#">MY100X-806</a> \$3,467.00	✓ 3-5 Days

Micropositioner

